

PE anti-human CD25 Antibody

| Catalog# / Size | 356103 / 25 tests 356104 / 100 tests 356134 / 100 μg |
|-------------------|---|
| Clone | M-A251 |
| Regulatory Status | RUO |
| Workshop | IV A053 |
| Other Names | IL-2 receptor α chain, Low affinity IL-2R, IL-2Rα chain |
| lsotype | Mouse IgG1, κ |
| Description | CD25 is a 55 kD type I transmembrane glycoprotein also known as low affinity IL-2 receptor α chain or Tac. It is expressed on progenitor lymphocytes, activated T and B cells, and activated monocytes/macrophages. CD25 is also expressed on a subset of non-stimulated CD4 ⁺ T cells termed T regulatory cells. Soluble CD25/IL-2R α is produced as a consequence of lymphocyte stimulation and is found in biological fluids following inflammatory responses. CD25 associates with IL-2 receptor β (CD122) and common γ (CD132) chains to form a high affinity IL-2R complex. |

Product Details

| Verified Reactivity | Human |
|---------------------|--|
| Reported Reactivity | Baboon, Cynomolgus, Rhesus |
| Antibody Type | Monoclonal |
| Host Species | Mouse |
| Immunogen | Human PHA-induced lymphocyte cells |
| Formulation | μg size: Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide. test sizes: Phosphate-buffered solution, pH 7.2, containing 0.09% sodium azide and BSA (origin USA). |
| Preparation | The antibody was purified by affinity chromatography and conjugated with PE under optimal conditions. |
| Concentration | μg sizes: 0.2 mg/mL test sizes: lot-specific (to obtain lot-specific concentration and expiration, please enter the lot number in our <u>Certificate of Analysis</u> online tool.) |
| Storage & Handling | The antibody solution should be stored undiluted between 2°C and 8°C, and protected from prolonged exposure to light. Do not freeze. |
| Application | FC - Quality tested |
| Recommended Usage | Each lot of this antibody is quality control tested by <u>immunofluorescent staining with flow</u> <u>cytometric analysis</u> . For flow cytometric staining using the μ g size, the suggested use of this reagent is $\leq 0.4 \ \mu$ g per million cells in 100 μ l volume. For flow cytometric staining using tests sizes, the suggested use of this reagent is 5 μ l per million cells in 100 μ l staining volume or 5 μ l per 100 μ l of whole blood. |
| Excitation Laser | Blue Laser (488 nm) Green Laser (532 nm)/Yellow-Green Laser (561 nm) |
| Application Notes | Additional reported applications (for the relevant formats) include: immunohistochemical staining of paraformaldehyde fixed frozen sections ¹ and spatial biology (IBEX) ^{2,3} . |
| | The CD25 molecule reveals three epitope regions: A, B, and C. M-A251 antibody recognizes |

epitope region B. Unlike other CD25 antibody clones, M-A251 can detect CD25 after fixation with paraformaldehyde.

| Application References (PubMed link indicates | 1. Li H and Pauza CD. 2015. <i>Eur. J. Immunol.</i> 45:298. (IHC) 2. Radtke AJ, <i>et al.</i> 2020. <i>Proc Natl Acad Sci USA</i> . 117:33455-33465. (SB) <u>PubMed</u> 3. Radtke AJ, <i>et al.</i> 2022. <i>Nat Protoc</i> . 17:378-401. (SB) PubMed |
|--|---|
| BioLegend citation) | 3. Native A3, et al. 2022. Nati 10:00. 11:310-401. (0.0) <u>rabited</u> |
| Product Citations | Chauvin C, <i>et al.</i> 2023. iScience. 26:106124. <u>PubMed</u> Cignarella F <i>et al.</i> 2018. Cell metabolism. 27(6):1222-1235. <u>PubMed</u> Parasar P, <i>et al.</i> 2022. Am J Reprod Immunol. 88:e13614. <u>PubMed</u> Ke F, <i>et al.</i> 2023. Elife. 12:. <u>PubMed</u> Deering RP, <i>et al.</i> 2023. Sci Rep. 13:8452. <u>PubMed</u> Serra-Peinado C, <i>et al.</i> 2019. Nat Commun. 10:3705. <u>PubMed</u> Schmitz T, <i>et al.</i> 2021. J Reprod Immunol. 148:103424. <u>PubMed</u> Beatson RE, <i>et al.</i> 2021. Cell Rep Med. 2:100473. <u>PubMed</u> Zhou Y, <i>et al.</i> 2017. Front Cell Infect Microbiol. 7:457. <u>PubMed</u> Holbrook AK, <i>et al.</i> 2019. Physiol Rep. 7:e14234. <u>PubMed</u> Ackerley CG, <i>et al.</i> 2022. Front Immunol. 13:972170. <u>PubMed</u> |
| RRID | AB_2561860 (BioLegend Cat. No. 356103) AB_2561861 (BioLegend Cat. No. 356104) AB_2564145 (BioLegend Cat. No. 356134) |

Antigen Details

| Structure | Type I transmembrane glycoprotein, 55 kD; low-affinity IL-2 receptor α chain |
|--------------------|---|
| Distribution | Activated T and B cells, monocytes/macrophages, Tregs |
| Interaction | Associates with IL-2R β /CD122 and IL-2R γ /CD132 receptor chains to form a high-affinity IL-2R complex |
| Ligand/Receptor | IL-2 |
| Cell Type | B cells, Macrophages, Monocytes, T cells, Tregs |
| Biology Area | Cell Biology, Immunology, Neuroscience, Neuroscience Cell Markers |
| Molecular Family | CD Molecules, Cytokine/Chemokine Receptors |
| Antigen References | Knapp W, <i>et al.</i> 1989. Leucocyte Typing IV: White Cell Differentiation Antigens. Oxford University Press. Schlossman S, <i>et al.</i> 1995. Leucocyte Typing V: White Cell Differentiation Antigens. Oxford University Press. Barclay N, <i>et al.</i> 1997. The Leukocyte Antigen FactsBook. Academic Press Inc. Taniguchi T and Minami Y. <i>et al.</i> 1993. <i>Cell</i> 73:5. Waldmann T. 1991. <i>J. Biol. Chem.</i> 266:2681. |
| Gene ID | 3559 |


Related Protocols

<u>Cell Surface Flow Cytometry Staining Protocol</u>

Other Formats

APC/Cyanine7 anti-human CD25, Purified anti-human CD25, PE anti-human CD25, FITC anti-human CD25, PE/Cyanine7 antihuman CD25, APC anti-human CD25, PerCP/Cyanine5.5 anti-human CD25, Brilliant Violet 421[™] anti-human CD25, Alexa Fluor® 488 anti-human CD25, Alexa Fluor® 700 anti-human CD25, Brilliant Violet 510[™] anti-human CD25, PE/Dazzle[™] 594 anti-human CD25, Biotin anti-human CD25, Alexa Fluor® 647 anti-human CD25, Pacific Blue[™] anti-human CD25, PerCP anti-human CD25, APC/Fire[™] 750 anti-human CD25, Brilliant Violet 711[™] anti-human CD25, Brilliant Violet 785[™] anti-human CD25, Brilliant Violet 605[™] anti-human CD25, KIRAVIA Blue 520[™] anti-human CD25, PE/Fire[™] 700 anti-human CD25, APC/Fire[™] 810 anti-human CD25, Spark NIR[™] 685 anti-human CD25 Antibody, Spark YG[™] 581 anti-human CD25, PE/Fire[™] 640 anti-human CD25 Antibody, PE/Cyanine5 anti-human CD25, Spark Red[™] 718 anti-human CD25, GMP PE anti-human CD25, GMP PE/Cyanine7 anti-human CD25, PerCP/Fire[™] 780 anti-human CD25, PE/Fire[™] 744 anti-human CD25, PerCP/Fire[™] 806 anti-human CD25, Brilliant Violet 650[™] anti-human CD25, Brilliant Violet 750[™] anti-human CD25

Product Data

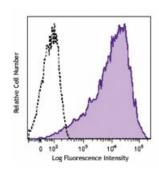

PHA-stimulated (3 day) human peripheral blood lymphocytes were stained with CD25 (clone M-A251) PE (filled histogram) or mouse IgG1, κ PE isotype control (open histogram).

For Research Use Only. Not for diagnostic or therapeutic use.

This product is supplied subject to the terms and conditions, including the limited license, located at www.biolegend.com/terms) ("Terms") and may be used only as provided in the Terms. Without limiting the foregoing, BioLegend products may not be used for any Commercial Purpose as defined in the Terms, resold in any form, used in manufacturing, or reverse engineered, sequenced, or otherwise studied or used to learn its design or composition without express written approval of BioLegend. Regardless of the information given in this document, user is solely responsible for determining any license requirements necessary for user's intended use and assumes all risk and liability arising from use of the product. BioLegend is not responsible for patent infringement or any other risks or liabilities whatsoever resulting from the use of its products.

BioLegend, the BioLegend logo, and all other trademarks are property of BioLegend, Inc. or their respective owners, and all rights are reserved.

8999 BioLegend Way, San Diego, CA 92121 www.biolegend.com Toll-Free Phone: 1-877-Bio-Legend (246-5343) Phone: (858) 768-5800 Fax: (877) 455-9587